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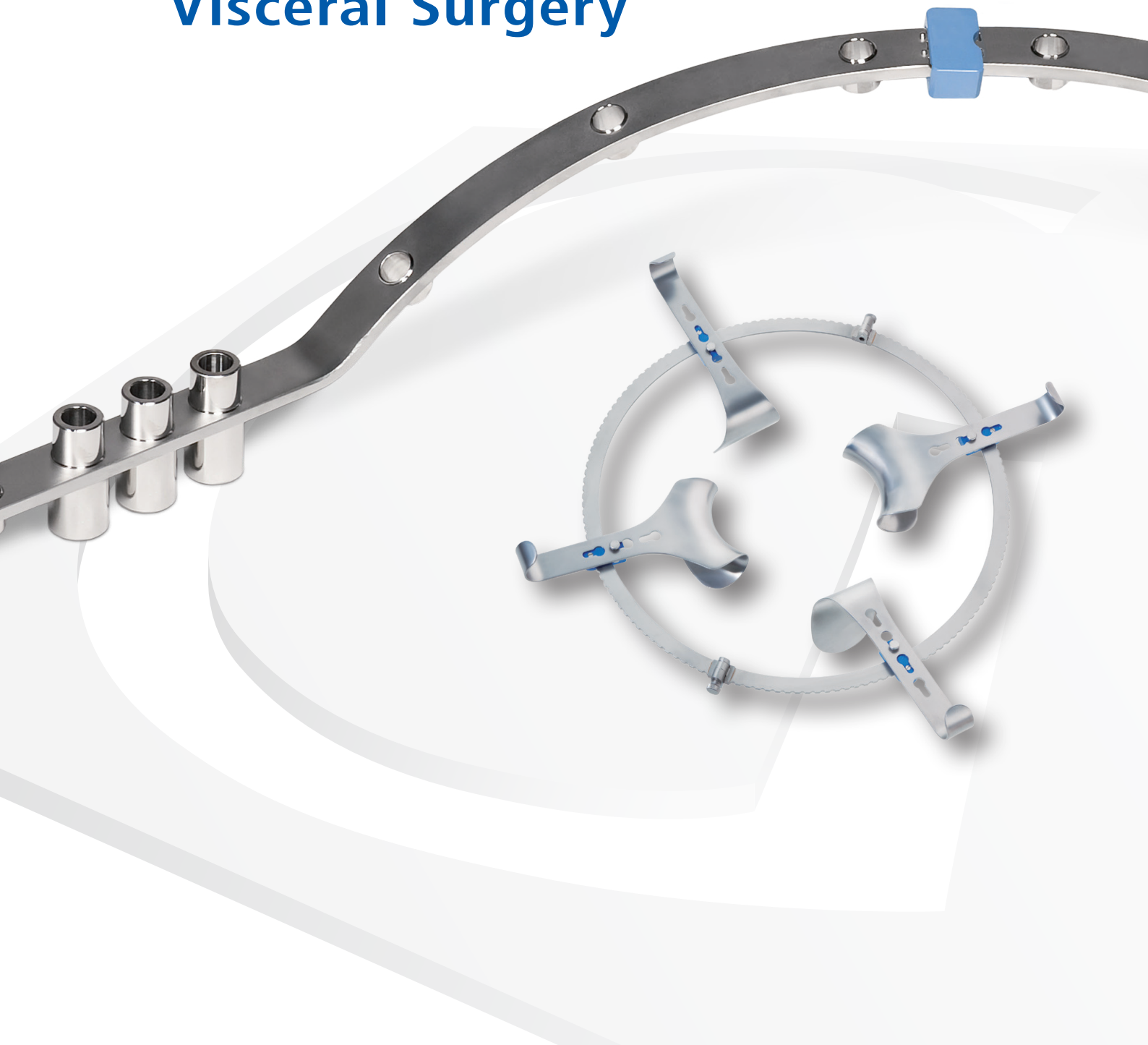
Our surgical instruments represent the successful result of many years of experience in medical technology – and quality made in Germany. We do our best every day to provide products that benefit our customers and patient health.

G BLUE LABEL™ is a comprehensive line of high-quality retractor systems for visceral surgery. In addition to conforming to a high standard of quality, G BLUE LABEL™ products are also easy to use and highly intuitive.

Rounding out the G BLUE LABEL™ are flexibility, fast delivery and personal assistance from our competent customer support staff.

For further information please contact

GEISTER Medizintechnik GmbH
Föhrenstrasse 2 | 78532 Tuttlingen/Germany
Tel: +49 74 61/966 240 | Fax: +49 74 61/966 24 22
www.gbluelabel.de

or visit the GEISTER® Extranet – Your portal for quality services at www.gbluelabel.de and download all catalogs as PDF files.

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If no other specifications are given all sizes mentioned in this catalog are in millimeter resp. inch.

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Infinitely variable retraction with the cable winch retractor

For precise, individual applications in the peritoneal cavity

Unique throughout the world, the cable winch technology used in the system allows surgeons to adjust retraction as desired with no preset positions when working in the peritoneal cavity. Thanks to flexible positioning for the insulated cables on the overhead bar, surgeons using the cable winch retractor can adapt the system to the situation at hand, adjusting the direction of traction and the position of the blades with exceptional precision. The retractor system is easy to mount on the upper section of the operating table, providing an optimum view of the operating site – even in difficult cases – and giving the surgeon the best possible visibility and range of motion.

NEW


Drehknopf aus PEEK
Adjusting knob made of PEEK

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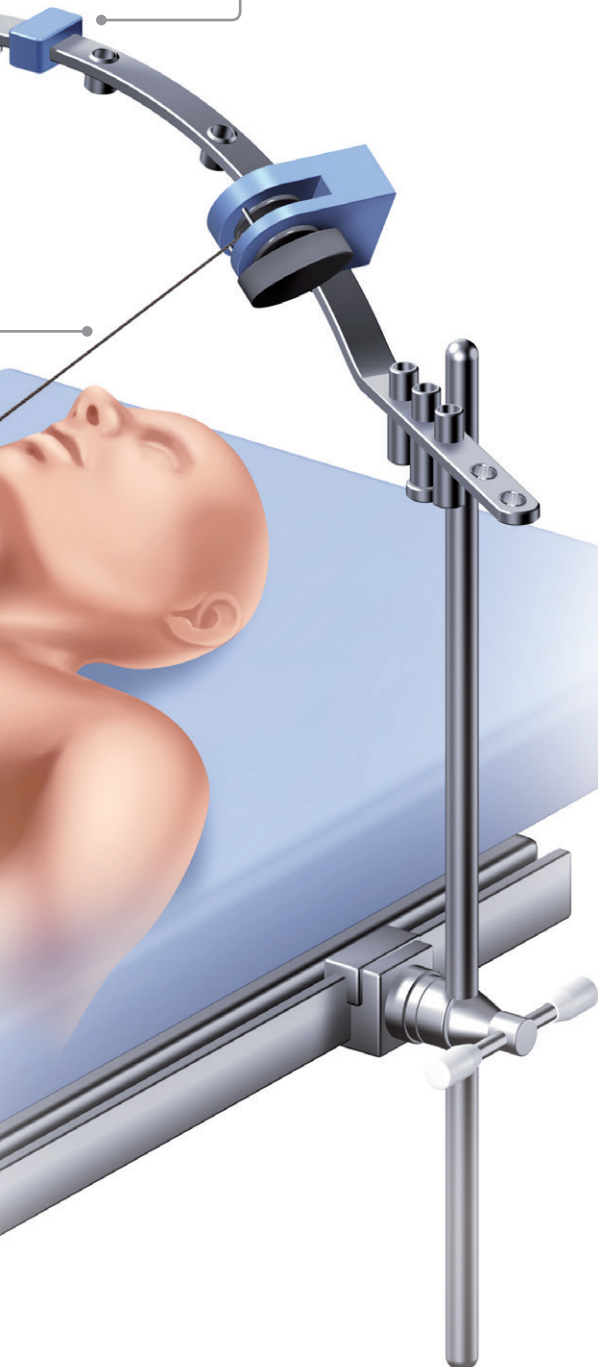- Ideal patient care and reduced tissue trauma thanks to a wide selection of blades in a variety of shapes and sizes
- The retractor is easy to set up and position thanks to the scalable support rods and the three mounting positions for the overhead bar – suitable for operating tables of any width
- The system can be processed in a standard sterilization container thanks to the foldable design of the frame
- Considerable stability and longevity thanks to carefully selected, high-quality materials
- Insulated cables
- NEW: Now features an extra-wide overhead bar for extreme lateral cable positioning, as well as larger blades for optimum care of obese patients.

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Innovative technical improvements for model from 2018

NEW

The cable is now completely coated

- Durable and reduces costs for cable exchange
- Higher tensile load of minimum 70 kg also for obese patients

The new material of the adjusting knob is PEEK

- Avoids breakage also in case of shock of frequent sterilization

The hinge at the frame is fixed by a PEEK slider. The pin is omitted.

- No lost pins
- Easy and fast mounting of the frame

All components of the Cable Winch Retractor are easy and safe in cleaning and processing. All improved components are compatible to existing systems.

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Benefits at a glance:

Optimized access to the operation field

- The lift of the thorax provides excellent view into the operation field
- The set up of the retractor is self-explanatory and provides an access for precise surgery
- Less assistance during surgery is needed

Simple and easy in handling


- Easy processing
- Simple set-up
- Infinitely variable retraction
- Traction can be applied from an extreme lateral position

Compact system

- Only essential components to perform a liver transplantation
- High durability of the system
- Further blades and frames available

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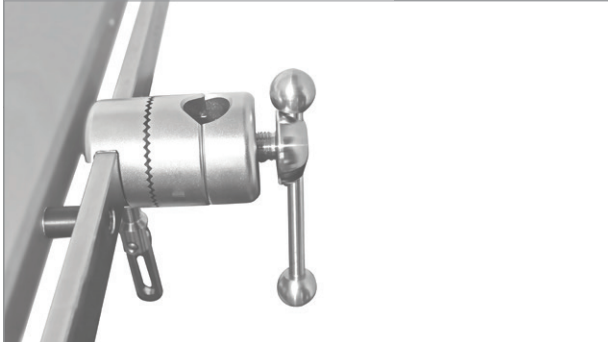1. The fixation element (UL7807-11) is secured to the standard rail of the OR table.

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2. The mounting support (UL7805-61/UL7805-62) is secured in the fixation element (UL7807-11).

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3. The offset blocks of the mounting support (UL7805-62) are mounted.

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4. The frame (UL7804-63 or UL7804-80) is placed on the mounting support.

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5. The cable winches (UL7801-09) are placed on the frame and provided with a blade.

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Both ends of the retractor system attach to the lateral bars of the operating table. The system consists of the following:

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**Cable winch Retractor
Complete set with blades**

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
**Cable winch Retractor extra wide
Complete set with blades**

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The system can be processed in a standard container and tray.

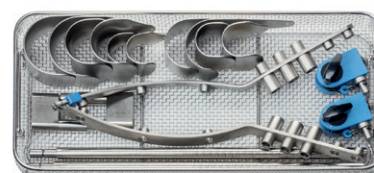GEISTER®

1

Seilzugsystem
Cable Winch System

Länge length	Breite width	Tiefe depth	Artikel-Nr. Product No.
88	63	69	UL7801-09

Feinstufig justierbar
Infinitely variable adjustment

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Reparatur-Set
Repair-Kit

Artikel-Nr. Product No.
UL7801-63

Inbusschlüssel und Loctite 5g
Hex-wrench and Loctite 5g

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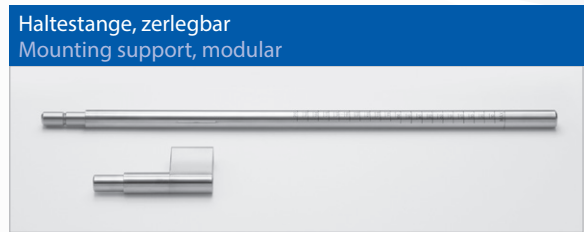
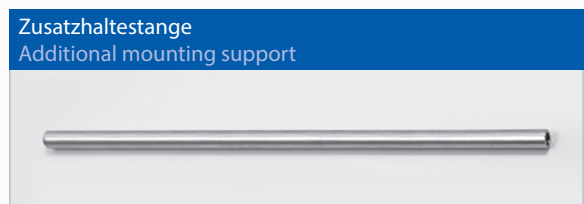
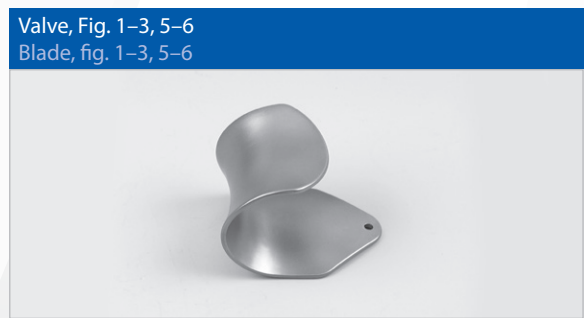
Classic visceral surgery with the rib retractor acc. to Rochard

For standard applications or surgeries of the upper abdomen

The rib retractor acc. to Rochard is considered to be the classic retractor for open stomach surgery. Because it applies constant traction, the system provides uninterrupted visibility of the operating site, up to the edge of the rib cage, providing ideal support for the user during an operation.

Because it is mounted on an insulated traction bar, the hook acc. to rochard can be adjusted to a variety of positions. The stand's two supports provide additional stability.

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
- Suitable for variable applications
thanks to a selection of five different hook widths
- A variety of positions
thanks to the hook acc. to rochard and traction bar, which combine to offer nine different positioning options
- Efficient processing
thanks to the removable horizontal rod and the resulting modest space requirements
- Considerable stability and longevity
thanks to carefully selected, high-quality materials
- Insulated traction bar

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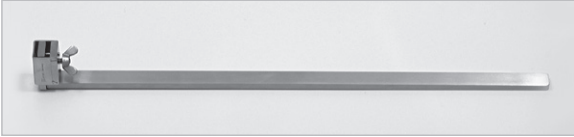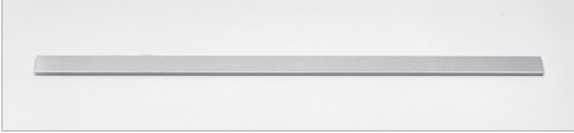
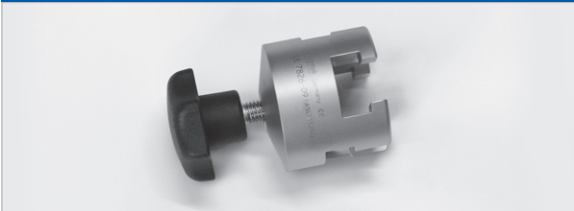
The established, efficient all-rounder – the rib retractor acc. to Behrends

For applications in the abdominal area

Because its design is simple and it is easy to use, the rib retractor acc. to Behrends has been used for years as a reliable partner in hospitals throughout the world.

Retractor hook positioning is flexible and secure, making the system suitable for use in a wide variety of applications, assisting the user even in difficult abdominal surgeries.

The rib retractor acc. to Stuhler may be added to the system.

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- Flexible use
thanks to a variety of retractors for different operating site sizes
- Various options for positioning retractor hooks and applying traction
thanks to mounting pins
- Can be processed in a standard sterilization container
- Considerable stability and longevity
thanks to carefully selected, high-quality materials
- Insulated attachment points for mounting bar

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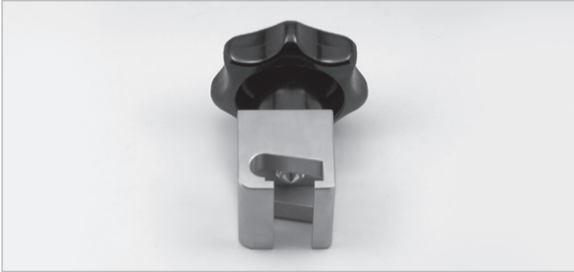
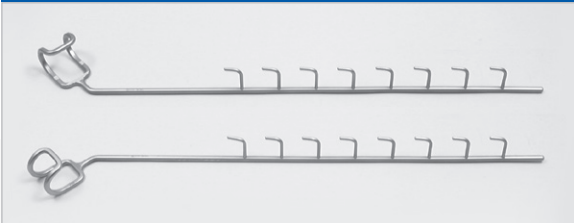
Simple and adaptable: The abdominal frame acc. to Zenker

For all kinds of abdominal surgery – suitable for use in combination with retractor systems

Whether used on its own or in combination with the retractor systems, the performance of the abdominal frame acc. to Zenker is impressive in a wide array of visceral surgery applications. Plus, this ingenious system has also been shown to be a helpful surgical tool in the fields of gynecology and urology.

The frame comes with four tabs for mounting and securely positioning blades through mounting holes in the blades. The tabs can be slid along the frame, making it easy to position the blades precisely, even when surgery is already underway. Blades with mounting hooks can also be used, and any number of these blades can be affixed to the frame at any location where no tabs are located.

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- The perfect frame size for every application thanks to a large selection of frame diameters – from 250 to 500 mm
- Can be processed in a standard sterilization container thanks to a foldable design
- Prevents additional tissue trauma thanks to a wide variety of blades

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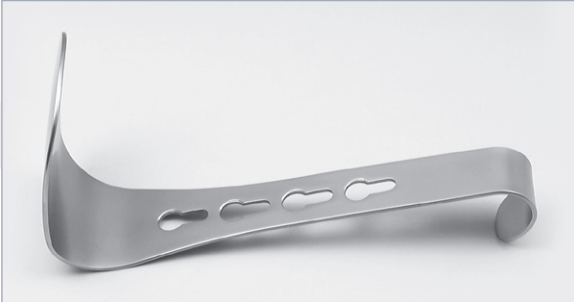
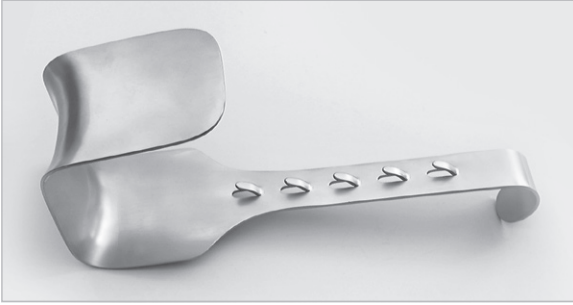
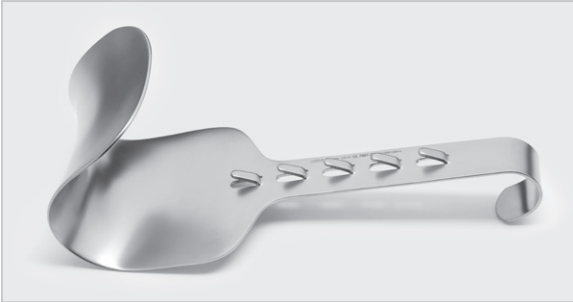
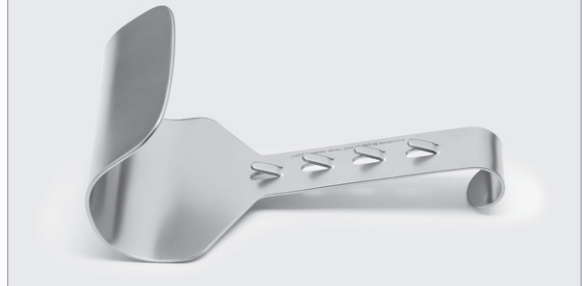
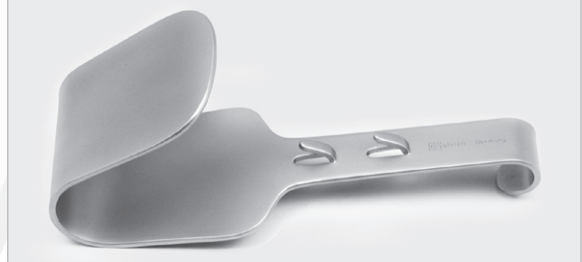
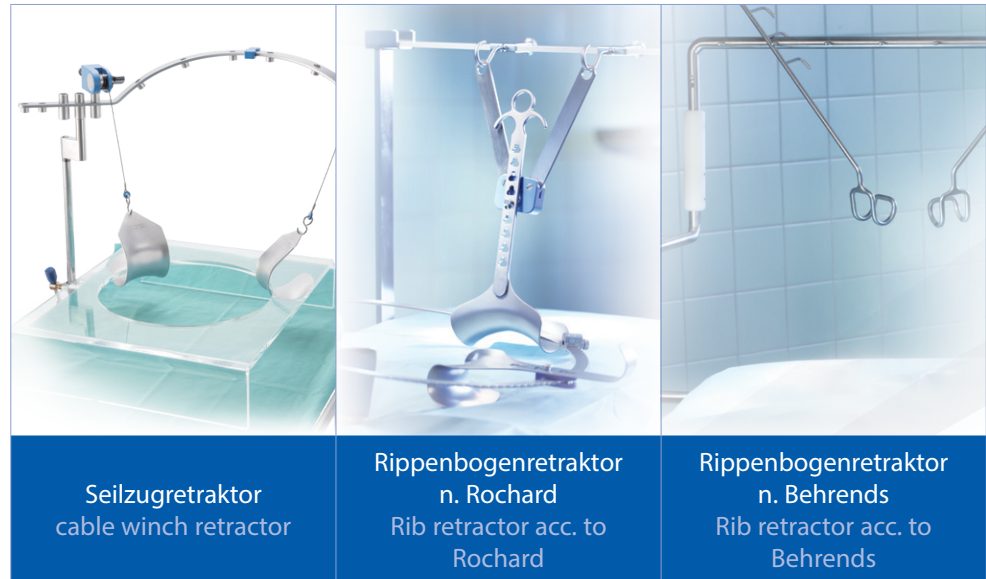
* when used with the rib retractor acc. to Stuhler


6

Artikelnummer Product number	Seite Page
UL7622-01-31	22
UL7622-01-36	22
UL7622-01-39	22
UL7622-01-44	22
UL7622-01-50	22
UL7622-25	22
UL7622-31	22
UL7622-36	22
UL7622-39	22
UL7622-44	22
UL7642-04	22
UL7661-20	22
UL7662-19	22
UL7663-21	23
UL7664-19	23
UL7665-22	23
UL7666-21	24
UL7667-22	24
UL7669-20	24
UL7670-20	23
UL7671-20	23
UL7671-26	23
UL7672-21	23
UL7673-19	24
UL7674-21	25
UL7675-18	24
UL7676-18	24
UL7677-19	25
UL7685-26	23
UL7690-20	25
UL7691-21	25
UL7695-17	25
UL7696-16	25
UL7800-00	9
UL7800-01	9
UL7801-09	9-10
UL7801-62	10
UL7801-63	10
UL7804-02	10
UL7804-63	9-10
UL7804-80	9-10
UL7805-61	9, 11
UL7805-62	9, 11

Artikelnummer Product number	Seite Page
UL7806-55	11
UL7807-11	9, 11
UL7810-39	9, 11
UL7811-42	9, 11
UL7812-55	9, 11
UL7813-80	9, 11
UL7814-80	11
UL7815-80	11
UL7821-67	14
UL7822-63	14
UL7823-12	14
UL7824-18	14
UL7825-68	14
UL7826-09	14
UL7828-10	15
UL7828-11	15
UL7828-12	15
UL7828-14	15
UL7828-16	15
UL7841-69	18
UL7842-11	18
UL7845-40	18
UL7845-50	18
UL7846-40	18
UL7846-50	18
UL7860-09	18
UL7860-12	18

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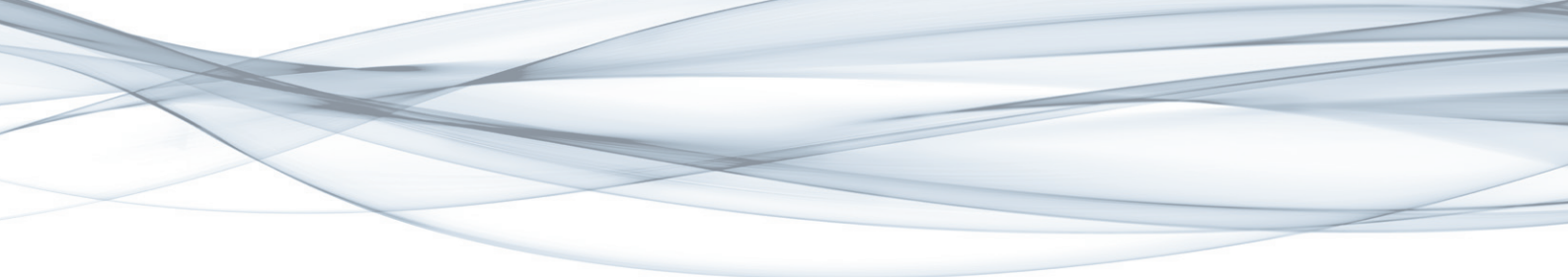Föhrenstrasse 2 | D-78532 Tuttlingen | Germany
☎ +49 7461/966240 | 📠 +49 7461/9662422
www.gbluelabel.de | info@gbluelabel.de