




## blunt fill needle

To be used with glass ampules with rubber stopper and plastic ampules

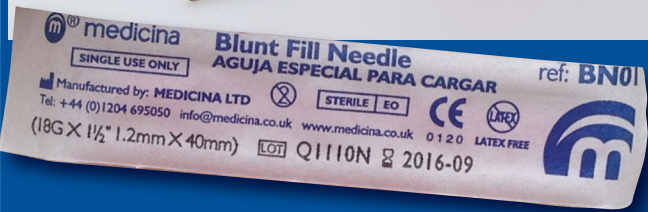

Gauge: 18G x 1.5" (38mm)  
 Packaging: Blue Writing  
 Hub: Red  
 Needles Shield: Red  
 Ref: **BN01**      NPC: **FTR1770**

## blunt filter needle

To be used with glass ampules only



Gauge: 18G x 1.5" (38mm)  
 Packaging: Green Writing  
 Hub: Lilac  
 Needles Shield: Red  
 Ref: **BN02**      NPC: **FTR1769**



When adding drugs to an infusion from a glass ampoule, please draw up using a **filter needle** and please remember to change to a **blunt needle** to inject into the bag.

Be sure to puncture the stopper centrally and at 90°. This will minimise the risk of catching the internal wall of the stopper.

*Innovative safety solutions*

